PROCISE/PROCISE cLC PROTEIN SEQUENCER

Site Preparation and Safety Guide

© Copyright 2002, Applied Biosystems

For Research Use Only. Not for use in diagnostic procedures.

ABI PRISM, the ABI PRISM design, Applied Biosystems, Aquapore, AmpliCover, Anitron, Biobytes, Brownlee, FastPhoramidite, GeneScan, Genotyper, HLP, INHERIT, Masterpiece, MicroAmp, MicroCoat, MPLC, NEWGUARD, ONESTEP, OPC, PCR-MATE, Phosphalink, POLYPORE, Precipitette, ProBlott, PROCISE, PROCISE cLC, ProFocus, ProSort, ProSpin, SeqEd, Sequence Navigator, SPHERIS, SPHERI10, StockMarks, Stretch, Synergy, SynthAssist, and VeloSep are registered trademarks of Applied Biosystems Corporation or its subsidiaries in the U.S. and certain other countries.

ABI, AmpFLSTR, AutoAssembler, BaseSprinter, CATALYST, GeneAssist, LV40, MatchMaker, PDQ, Primer Express, and ProSorb are trademarks of Applied Biosystems Corporation or its subsidiaries in the U.S. and certain other countries.

All other trademarks are the sole property of their respective owners.

Contents

1 Site Preparation

Before You Start.....	1-1
Preinstallation	1-1
Operator Training.....	1-1
Performance Verification	1-1
Ordering Supplies	1-1
Unpacking	1-1
Technical Support.....	1-2
To Reach Us On the Web.....	1-2
To Contact Technical Support by E-Mail	1-2
Hours for Telephone Technical Support	1-2
To Contact Technical Support by Telephone or Fax	1-2
To Reach Technical Support Through the Internet	1-5
To Obtain Documents-on-Demand	1-5
Preinstallation Checklist.....	1-6
About this Checklist	1-6
Checklist	1-6
Items Shipped with this Instrument	1-8
Summary	1-8
Hazardous Chemicals Warning	1-8
Installation Chemical Kit.....	1-9
Items Needed but Not Supplied	1-10
Resupply	1-10
Water	1-10
Equipment and Supplies	1-10
Laboratory Safety	1-11
Onsite Representative	1-11
Required Safety Equipment.....	1-11
Laboratory Space Required	1-12
Dimensions	1-12
Location of Instrument	1-12
Clearance	1-12
Typical Laboratory Layout	1-13
Electrical Requirements	1-14
Power	1-14
Grounding	1-14
Power Cords.....	1-14
Voltage Quality	1-14

Voltage Spikes	1-14
Power Outages	1-15
Laboratory Environmental Requirements	1-16
Altitude	1-16
Temperature and Humidity	1-16
Pollution	1-16
Heat	1-16
Emission/ Immunity Statement	1-16
Laboratory Ventilation	1-17
Venting	1-17
Fume Hood	1-17
Duct System	1-17

2 Instrument Safety

Instrument Safety	2-1
Safe Operation	2-1
User Attention Words	2-1
Instrument Labeling	2-1
Safety Labels	2-1
Signal Words	2-1
Labels That May Be Found On Your Instrument	2-2
Safety Alert Symbols	2-3
Electrical Symbols	2-3
Non-electrical Symbols	2-3
Input/Output Connections	2-4
Location	2-4
Gas Safety	2-5
Pressurized Gas Safety	2-5
Gaseous Waste Safety	2-5
Disposing of Gaseous Waste	2-6
Instrument Waste and Exhaust System	2-6
Exhaust Line Tubing	2-6
Connecting the Exhaust line	2-6
Pressurized Gas and Accessories Needed but Not Supplied	2-7
Gas Cylinder	2-7
Pressure Regulator	2-7
Attaching the Cylinder	2-7

3 Chemical Safety

Introduction	3-1
Please read...	3-1

Hazardous Chemicals	3-2
Overview	3-2
Handling Hazardous Chemicals	3-2
Hazardous Waste	3-3
Overview	3-3
Handling Hazardous Waste	3-3
Storing Hazardous Waste	3-3
Instrument Waste Profile Overview	3-4
General Information.	3-4
Waste Profile Sections	3-4
Material Safety Data Sheets	3-9
Overview	3-9
Contents of MSDSs	3-9
How to Get Extra Copies of MSDSs	3-9
MSDSs for the PROCISE cLC Protein Sequencer.	3-10
MSDSs for Applied Biosystems Chemicals	3-10
MSDSs for Chemicals Not Supplied by Applied Biosystems	3-10

A Acronyms and Abbreviations

Acronyms and Abbreviations Used in Waste Profiles and MSDSs	A-1
Introduction	A-1
Organizations, Regulations, and Scientific Terminology	A-1
Units of Measurement	A-2
Chemicals	A-3

Site Preparation

1

Before You Start...

Preinstallation	Before the instrument is installed, the installation site must be prepared so that the instrument can operate correctly and safely. Careful attention to the requirements presented here will simplify the installation procedure.
------------------------	--

Operator Training	Training of operators is one of the primary goals of installation. Persons who are to be trained to operate the instrument should set aside two uninterrupted days to work with the Applied Biosystems service representative. If this is not possible, the installation should be rescheduled.
--------------------------	---

Performance Verification	Calibration of the instrument and verification of performance will be performed by an Applied Biosystems service representative during the installation.
---------------------------------	--

Ordering Supplies	Reagents shipped with this instrument will be consumed during the process of setup and verification. Before installation, be sure to order additional chemicals and other supplies necessary for the ongoing operation of the instrument.
--------------------------	---

Unpacking	Do not unpack instruments. Inspect instrument cartons and report any damage to your Applied Biosystems service representative. For information and instructions about unpacking the Installation Chemical Kit, see "Items Shipped with this Product" later in this chapter.
------------------	---

Technical Support

Contacting Technical Support

You can contact Applied Biosystems for technical support by telephone or fax, by e-mail, or through the Internet. You can order Applied Biosystems user documents, MSDSs, certificates of analysis, and other related documents 24 hours a day. In addition, you can download documents in PDF format from the Applied Biosystems Web site (please see the section "To Obtain Documents on Demand" following the telephone information below).

To Contact Technical Support by E-Mail

Contact technical support by e-mail for help in the following product areas:

Product Area	E-mail address
Genetic Analysis (DNA Sequencing)	galab@appliedbiosystems.com
Sequence Detection Systems and PCR	pcrlab@appliedbiosystems.com
Protein Sequencing, Peptide and DNA Synthesis	corelab@appliedbiosystems.com
Biochromatography, PerSeptive DNA, PNA and Peptide Synthesis systems, CytoFluor [®] , FMAT [™] , Voyager [™] , and Mariner [™] Mass Spectrometers	tsupport@appliedbiosystems.com
LC/MS (Applied Biosystems/MDS Sciex)	apisupport@sciex.com or api3-support@sciex.com
Chemiluminescence (Tropix)	tropix@appliedbiosystems.com

Hours for Telephone Technical Support

In the United States and Canada, technical support is available at the following times:

Product	Hours
Chemiluminescence	8:30 a.m. to 5:30 p.m. Eastern Time
Framingham support	8:00 a.m. to 6:00 p.m. Eastern Time
All Other Products	5:30 a.m. to 5:00 p.m. Pacific Time

To Contact Technical Support by Telephone or Fax

In North America

To contact Applied Biosystems Technical Support, use the telephone or fax numbers given below. (To open a service call for other support needs, or in case of an emergency, dial **1-800-831-6844** and press **1**.)

Product or Product Area	Telephone Dial...	Fax Dial...
ABI PRISM [®] 3700 DNA Analyzer	1-800-831-6844, then press 8	1-650-638-5981
DNA Synthesis	1-800-831-6844, then press 21	1-650-638-5981
Fluorescent DNA Sequencing	1-800-831-6844, then press 22	1-650-638-5981
Fluorescent Fragment Analysis (includes GeneScan [®] applications)	1-800-831-6844, then press 23	1-650-638-5981

Product or Product Area	Telephone Dial...	Fax Dial...
Integrated Thermal Cyclers (ABI PRISM® 877 and Catalyst 800 instruments)	1-800-831-6844, then press 24	1-650-638-5981
ABI PRISM® 3100 Genetic Analyzer	1-800-831-6844, then press 26	1-650-638-5981
BioInformatics (includes BioLIMS®, BioMerge™, and SQL GT™ applications)	1-800-831-6844, then press 25	1-505-982-7690
Peptide Synthesis (433 and 43X Systems)	1-800-831-6844, then press 31	1-650-638-5981
Protein Sequencing (Procise® Protein Sequencing Systems)	1-800-831-6844, then press 32	1-650-638-5981
PCR and Sequence Detection	1-800-762-4001, then press 1 for PCR, 2 for the 7700 or 5700, 6 for the 6700 or dial 1-800-831-6844, then press 5	1-240-453-4613
Voyager™ MALDI-TOF Biospectrometry and Mariner™ ESI-TOF Mass Spectrometry Workstations	1-800-899-5858, then press 13	1-508-383-7855
Biochromatography (BioCAD® Workstations and Poros® Perfusion Chromatography Products)	1-800-899-5858, then press 14	1-508-383-7855
Expedite™ Nucleic acid Synthesis Systems	1-800-899-5858, then press 15	1-508-383-7855
Peptide Synthesis (Pioneer™ and 9050 Plus Peptide Synthesizers)	1-800-899-5858, then press 15	1-508-383-7855
PNA Custom and Synthesis	1-800-899-5858, then press 15	1-508-383-7855
FMAT™ 8100 HTS System and Cytofluor® 4000 Fluorescence Plate Reader	1-800-899-5858, then press 16	1-508-383-7855
Chemiluminescence (Tropix)	1-800-542-2369 (U.S. only), or 1-781-271-0045	1-781-275-8581
Applied Biosystems/MDS Sciex	1-800-952-4716	1-650-638-6223

Outside North America

Region	Telephone Dial...	Fax Dial...
Africa and the Middle East		
Africa (English Speaking) and West Asia (Fairlands, South Africa)	27 11 478 0411	27 11 478 0349
South Africa (Johannesburg)	27 11 478 0411	27 11 478 0349
Middle Eastern Countries and North Africa (Monza, Italia)	39 (0)39 8389 481	39 (0)39 8389 493

Region	Telephone Dial...	Fax Dial...
Eastern Asia, China, Oceania		
Australia (Scoresby, Victoria)	61 3 9730 8600	61 3 9730 8799
China (Beijing)	86 10 64106608	86 10 64106617
Hong Kong	852 2756 6928	852 2756 6968
Korea (Seoul)	82 2 593 6470/6471	82 2 593 6472
Malaysia (Petaling Jaya)	60 3 758 8268	60 3 754 9043
Singapore	65 896 2168	65 896 2147
Taiwan (Taipei Hsien)	886 2 22358 2838	886 2 2358 2839
Thailand (Bangkok)	66 2 719 6405	66 2 319 9788
Europe		
Austria (Wien)	43 (0)1 867 35 75 0	43 (0)1 867 35 75 11
Belgium	32 (0)2 712 5555	32 (0)2 712 5516
Czech Republic and Slovakia (Praha)	420 2 61 222 164	420 2 61 222 168
Denmark (Naerum)	45 45 58 60 00	45 45 58 60 01
Finland (Espoo)	358 (0)9 251 24 250	358 (0)9 251 24 243
France (Paris)	33 (0)1 69 59 85 85	33 (0)1 69 59 85 00
Germany (Weiterstadt)	49 (0) 6150 101 0	49 (0) 6150 101 101
Hungary (Budapest)	36 (0)1 270 8398	36 (0)1 270 8288
Italy (Milano)	39 (0)39 83891	39 (0)39 838 9492
Norway (Oslo)	47 23 12 06 05	47 23 12 05 75
Poland, Lithuania, Latvia, and Estonia (Warszawa)	48 (22) 866 40 10	48 (22) 866 40 20
Portugal (Lisboa)	351 (0)22 605 33 14	351 (0)22 605 33 15
Russia (Moskva)	7 095 935 8888	7 095 564 8787
South East Europe (Zagreb, Croatia)	385 1 34 91 927	385 1 34 91 840
Spain (Tres Cantos)	34 (0)91 806 1210	34 (0)91 806 1206
Sweden (Stockholm)	46 (0)8 619 4400	46 (0)8 619 4401
Switzerland (Rotkreuz)	41 (0)41 799 7777	41 (0)41 790 0676
The Netherlands (Nieuwerkerk a/d IJssel)	31 (0)180 331400	31 (0)180 331409
United Kingdom (Warrington, Cheshire)	44 (0)1925 825650	44 (0)1925 282502
All other countries not listed (Warrington, UK)	44 (0)1925 282481	44 (0)1925 282509
Japan		
Japan (Hacchobori, Chuo-Ku, Tokyo)	81 3 5566 6230	81 3 5566 6507
Latin America		
Del.A. Obregon, Mexico	305-670-4350	305-670-4349

To Reach Technical Support Through the Internet

We strongly encourage you to visit our Web site for answers to frequently asked questions and for more information about our products. You can also order technical documents or an index of available documents and have them faxed or e-mailed to you through our site. The Applied Biosystems Web site address is

<http://www.appliedbiosystems.com/techsupp>

To submit technical questions from North America or Europe:

Step	Action
1	Access the Applied Biosystems Technical Support Web site.
2	Under the Troubleshooting heading, click Support Request Forms, then select the relevant support region for the product area of interest.
3	Enter the requested information and your question in the displayed form, then click Ask Us RIGHT NOW (blue button with yellow text).
4	Enter the required information in the next form (if you have not already done so), then click Ask Us RIGHT NOW. You will receive an e-mail reply to your question from one of our technical experts within 24 to 48 hours.

To Obtain Documents on Demand

Free, 24-hour access to Applied Biosystems technical documents, including MSDSs, is available by fax or e-mail or by download from our Web site.

To order documents...	Then...
by index number	<ol style="list-style-type: none"> Access the Applied Biosystems Technical Support Web site at http://www.appliedbiosystems.com/techsupp Click the Index link for the document type you want, then find the document you want and record the index number. Use the index number when requesting documents following the procedures below.
by phone for fax delivery	<ol style="list-style-type: none"> From the U.S. or Canada, call 1-800-487-6809, or from outside the U.S. and Canada, call 1-858-712-0317. Follow the voice instructions to order the documents you want. <p>Note There is a limit of five documents per request.</p>
through the Internet for fax or e-mail delivery	<ol style="list-style-type: none"> Access the Applied Biosystems Technical Support Web site at http://www.appliedbiosystems.com/techsupp Under Resource Libraries, click the type of document you want. Enter or select the requested information in the displayed form, then click Search. In the displayed search results, select a check box for the method of delivery for each document that matches your criteria, then click Deliver Selected Documents Now (or click the PDF icon for the document to download it immediately). Fill in the information form (if you have not previously done so), then click Deliver Selected Documents Now to submit your order. <p>Note There is a limit of five documents per request for fax delivery but no limit on the number of documents you can order for e-mail delivery.</p>

Preinstallation Checklist

About this Checklist Use this checklist to ensure that all preparations have been made for installation. A service representative will contact you to verify that everything is checked off before the installation date.

Checklist Check off and enter date to ensure that all preparations have been made.

PROCISE/PROCISE cLC Protein Sequencer Preinstallation Checklist

✓ if ready	Date Confirmed	Components
		General
		Received instrument(s) and verified that cartons are intact.
		Read this <i>Site Preparation and Safety Manual</i> , including MSDSs.
		Verified that instrument(s), serial number(s), and system configuration, as shown on the packing list, are the same as ordered.
		Set aside two uninterrupted days for in-lab training during installation.
		Unpacked and stored contents of Installation Chemical Kit.
		Electrical
		A dedicated 2.0 kVA power line and ground, or a 2.0 kVA power line with a line conditioner or UPS, is in place.
		At least four standard power outlets are within 2.5 m (8 ft.) of the instrument location.
		Instrument voltage, if specified on the packing list, matches the voltage available in the laboratory.
		Laboratory
		Laboratory bench is of correct dimensions to accommodate the system and is situated so that the instrument is accessible to the installer on all four sides.
		Laboratory safety requirements, as specified in this manual, have been met.
		Laboratory environmental requirements, as specified in this manual, have been met.
		Room ventilation accommodates instrument heat output.
		Deionized water is on site.
		Proper waste disposal method for hazardous chemical waste has been established.
		Fume hood or ventilation duct system is located within 5 m (15 ft.) of the instrument.
		Equipment
		Waste Bottle with secondary containment
		Macintosh-compatible printer; this cannot be ordered with the instrument
		A two-stage pressure regulator with dual gauges (output range: 0–80 psi, high pressure range: 0–3000 psi), and a Compressed Gas Association (CGA) 580 cylinder adapter for each gas cylinder

PROCISE/PROCISE cLC Protein Sequencer Preinstallation Checklist *(continued)*

✓ if ready	Date Confirmed	Components
		Approved strap or clamp to safely secure gas cylinders to the wall or bench/table
		Consumable Supplies
		Additional Applied Biosystems reagent kits
		Additional computer supplies (paper, disks, etc.)
		Chemically resistant, powder-free disposable gloves
		Two Argon gas cylinders of 99.998% purity or greater
		Acetone (HPLC-grade)
		Methanol (in wash bottle)
		Sodium Dodecyl Sulfate (SDS)

Items Shipped with this Instrument

Summary The PROCISE/PROCISE cLC Protein Sequencer is shipped with the following:

- ◆ User's manual
- ◆ Installation Chemical Kit (P/N 401690 and 402064)
- ◆ Apple Macintosh Power PC computer system

Do not unpack any cartons related to the instrument itself. This will protect you from liability if any damage occurred during shipping.

IMPORTANT You must unpack the Installation Chemical Kit and store the components as specified in Table 1-1.

Hazardous Chemicals Warning

! WARNING ! Some chemicals used with Applied Biosystems instruments are hazardous and can cause injury, illness, or death. Always read the appropriate Waste Profiles and MSDSs before interacting with the instrument and chemicals in any way. Hazardous Chemical Warnings are prominently displayed on the labels of all hazardous materials.

Before unpacking the Installation Chemical Kit, storing chemicals, or interacting with the chemicals and instrument, read Chapter 3, "Chemical Safety," of this manual. Chapter 3 contains the Waste Profiles and Material Safety Data Sheets (MSDSs) that pertain to this instrument.

continued on next page

Installation Chemical Kit Unpack the Installation Chemical Kit upon receipt. The Installation Chemical Kit is intended to be used during installation to verify instrument performance.

Store the chemicals and reagents as indicated in the table below.

Note All reagents are common to both kits below (P/N 401690 and P/N 402064) unless otherwise indicated.

Table 1-1 PROCISE/PROCISE cLC Protein Sequencer Installation Chemical Kit Components (P/N 401690 and P/N 402064)

P/N	Description	Storage Conditions
400208	R1, 5% PITC in N-Heptane	20 °C
401535	R2B, n-Methylpiperidine in methanol and water	4 °C
400003	R3, Trifluoroacetic Acid	Room temp
400028	R4A, 25% Trifluoroacetic Acid in water	4 °C
400315	R5, Acetonitrile with DTT	Room temp
400079	S1, N-Heptane (Procise cLC only)	Room temp
400854	S2B, Ethyl Acetate	Room temp
400008	S3, N-Butyl Chloride	Room temp
400314	S4B, 20% Acetonitrile in water (Procise only)	Room temp
400470	MeOH (Methanol)	Room temp
400879	PTH Analyzer Standards	20 °C
400385	Biobrene Plus	4 °C
400979	β-Lactoglobulin	4 °C
401446	Premix buffer concentrate	4 °C
402051	S4C, 10% Acetonitrile in water, 1 L (Procise)	Room temp
401464	Solvent A3, 3.5% THF Tetrahydrofuran in Water, 1 L (Procise)	Room temp
401570	Solvent B2, 12% Isopropanol and Acetonitrile, 1 L (Procise)	Room temp
401887	Solvent A3, 3.5% THF Tetrahydrofuran in Water, 450 mL (cLC)	Room temp
401886	Solvent B2, 12% Isopropanol and Acetonitrile, 450 mL (cLC)	Room temp

Note For the Procise sequencer, solvents will be in 1 L bottles. For the cLC sequencer, solvents will be in 450 mL bottles.

Items Needed but Not Supplied

Resupply	The Installation Chemical Kit components are completely consumed during the installation and initial testing of the instrument. To ensure an uninterrupted supply of reagents and eliminate the higher transportation costs of rush shipments, you should order additional chemicals and supplies in advance.
Water	Although the Installation Chemical Kit contains all the reagents necessary for operating the instrument, each laboratory must provide its own source of high-quality deionized water.
Equipment and Supplies	<p>In addition to the Installation Chemical Kit, the following equipment and supplies are needed for operation of the PROCISE/PROCISE cLC Protein Sequencer:</p> <ul style="list-style-type: none">◆ Waste bottle with secondary containment◆ Macintosh-compatible printer; this cannot be ordered with the instrument◆ A two-stage pressure regulator with dual gauges (output range: 0–80 psi; high pressure range: 0–3000 psi), and a Compressed Gas Association (CGA) 580 cylinder adapter for each gas cylinder◆ Approved strap or clamp to safely secure gas cylinders to the wall, bench, table, or cart◆ Additional Applied Biosystems reagent kits◆ Additional computer supplies (paper, disks, etc.)◆ Chemically resistant disposable gloves, lab coats, and safety glasses or goggles◆ Two Argon gas cylinders of 99.998% purity or greater◆ Acetone (HPLC-grade)◆ Methanol (in wash bottle)◆ Sodium Dodecyl Sulfate (SDS)

Laboratory Safety

Onsite Representative	We request that a representative from your laboratory be in the vicinity and available to the Applied Biosystems engineer at all times while he or she is on-site.
------------------------------	--

Required Safety Equipment	Your laboratory has specific safety practices and policies designed to protect laboratory personnel from potential hazards, both obvious and hidden, that are present. We expect that all applicable, safety-related procedures will be followed at all times.
----------------------------------	--

The following safety equipment should be available:

- ◆ fire extinguisher (Halon)
 - ◆ eye wash
 - ◆ safety shower
 - ◆ eye and hand protection
 - ◆ adequate ventilation
 - ◆ first aid equipment
 - ◆ spill clean-up equipment
 - ◆ protection from any sources of radiation (lasers, radioisotopes, contaminated equipment, radioactive wastes, etc.) that may be present in the area where our engineer will be working
-

Laboratory Space Required

Dimensions The PROCISE/PROCISE cLC Protein Sequencer and the Macintosh computer have the following dimensions:

Component	Width	Depth	Height	Total Weight
Instrument	126.0 cm (49 in.)	103.5 cm (41 in.)	111.5 cm (44 in.)	125 kg (275 lb.)
Computer	35.5 cm (14 in.)	43.0 cm (17 in.)	45.5 cm (18 in.)	

Location of Instrument The PROCISE/PROCISE cLC Protein Sequencer must be located within 5 m (15 ft.) of a fume hood or ventilation duct system. The instrument requires a laboratory bench or sturdy table of correct dimensions to accommodate the system. The instrument must be accessible to the installer on all four sides.

The Apple Macintosh computer and printer should be on the same table as the PROCISE/PROCISE cLC Protein Sequencer. If this is impractical, the computer must be located within 2 m (6 ft.) of the instrument's left side panel.

IMPORTANT The operator should be able to reach the instrument while viewing the computer screen.

Clearance A clearance of 30 cm (1 ft.) is needed at the rear of the instrument to provide adequate ventilation. Additional clearance is required for servicing. Do not block access to the rear of the instrument.

Sufficient space within approximately 1.5 m (5 ft.) of the instrument is needed to safely secure at least one size 1A Argon gas cylinder.

Note Additional information related to compressed gas safety issues is provided in Chapter 2, "Instrument Safety."

Sufficient space to the right side of the instrument is needed for solvent and waste bottles.

continued on next page

Typical Laboratory Layout A typical laboratory layout for the PROCISE/PROCISE cLC Protein Sequencer and accessory equipment is shown in Figure 1-1.

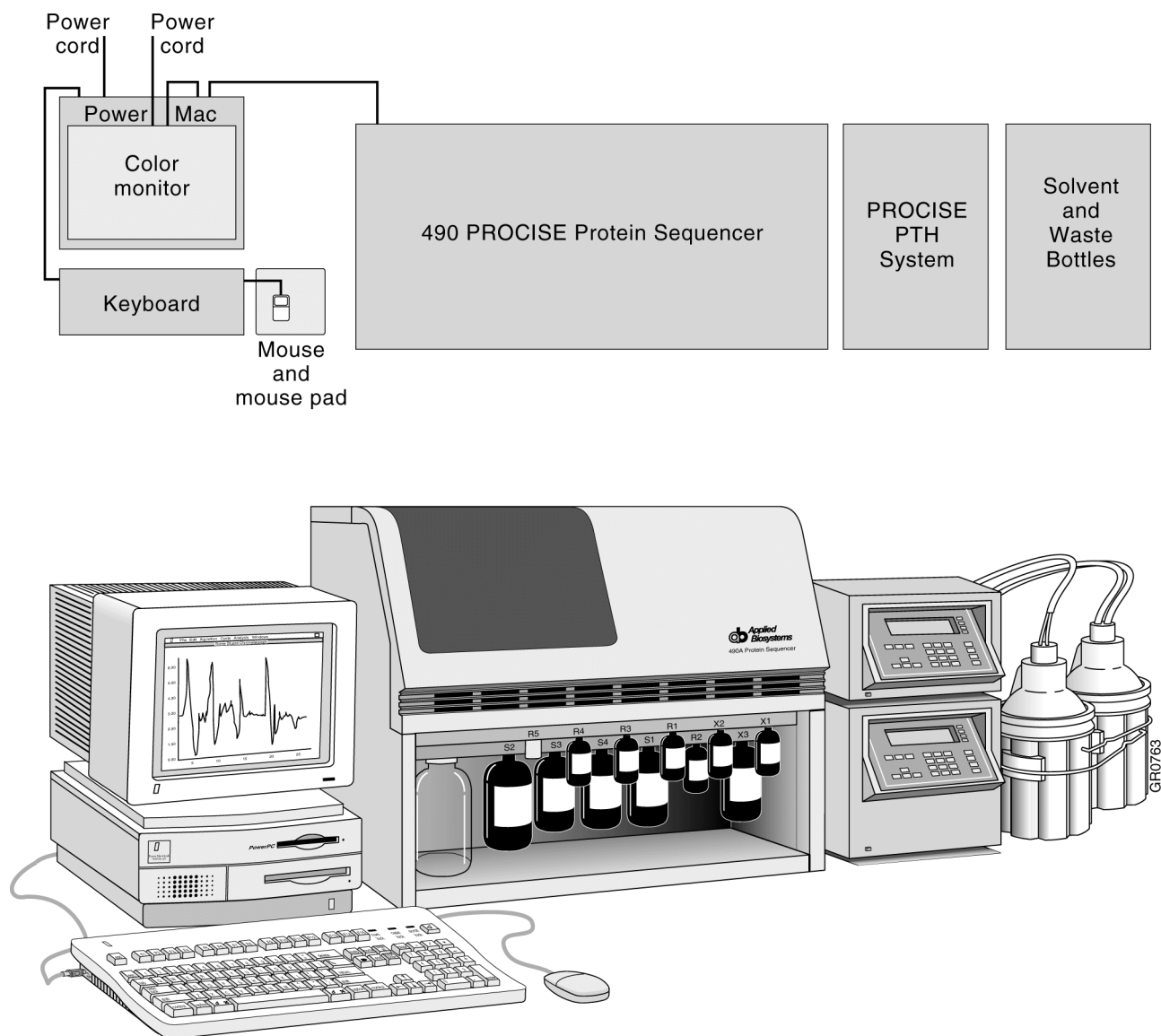
Figure 1-1 Typical laboratory layout for the PROCISE/PROCISE cLC Protein Sequencer Sequencing System

Electrical Requirements

Power The electrical receptacle should have a dedicated 2.0 kVA power line and ground or a 2.0 kVA power line with a line conditioner or uninterruptible power supply (UPS). The electrical receptacle must be located within 2.5 m (8 ft.) of the instrument rear panel. The following table specifies the electrical operating range for various parts of the world.

Location	Volts (AC)	Frequency
Japan	100 \pm 10%	50/60 Hz \pm 1%
USA/Canada	120 \pm 10%	50/60 Hz \pm 1%
Europe (pre-1992)	220 \pm 10%	50/60 Hz \pm 1%
EC	230 \pm 10%	50/60 Hz \pm 1%
UK (pre-1992)	240 + 6%/-10%	50/60 Hz \pm 1%
Australia	240 + 6%/-10%	50/60 Hz \pm 1%

Grounding Certain types of electrical noise are greatly exaggerated by poor or improper electrical ground connections. To prevent these problems, it is very important to have a dedicated line and ground between the instrument and building main electrical service.

Power Cords In the USA, Canada, and Japan, the instrument is supplied with a detachable cord equipped with a standard three-prong plug.

In Europe and Australia, the instrument is supplied with an detachable electrical cord equipped with a standard EC plug.

The Macintosh computer can be plugged into any standard electrical receptacle after it has been configured for the proper voltage.

Voltage Quality Line voltage must be within $\pm 10\%$ of the nominal value. High or low voltages may have adverse effects on the electronic components of the instrument. In areas where the supplied power is subject to fluctuations exceeding these limits, a power line regulator may be required.

Voltage Spikes Short-duration, high-voltage spikes often cause random failures in microprocessor controlled instrumentation. These spikes can be caused by other devices using the same power source (refrigerators, air conditioners, and centrifuges) or by outside influences such as lightning. A dedicated line and ground between the instrument and building main electrical service are necessary to prevent such problems.

If your environment contains devices that are electrically noisy or you are in an area with frequent electrical storms, a line conditioner with a recommended capacity of 2.0 kVA will enhance the reliability of your system. This may be lower depending on the conditioner or power supply design.

continued on next page

Power Outages The instrument has been designed to recover from short periods of power outage (loss) and continue operation, provided that the line voltage does not become excessively noisy before the outage. If increased protection during a power outage is desired, you may want to install an uninterruptible power supply (UPS). We recommend a capacity of 2.0 kVA. The UPS will involve a higher cost than a line conditioner.

Laboratory Environmental Requirements

Altitude	This instrument is for indoor use only and for altitudes not exceeding 2,000 meters (6,500 ft.) above sea level.
Temperature and Humidity	The laboratory temperature should be maintained between 15–30 °C (59–85 °F). The instrument can tolerate up to 80% relative humidity. Avoid placing the instrument adjacent to heaters or cooling ducts.
Pollution	The installation category (overvoltage category) for this instrument is II, and it is classified as portable equipment. The instrument has a pollution degree rating of 2, and may be installed in an environment that has non-conductive pollutants only.
Heat	The thermal output of the instrument is 8800 Btu/h (~2,600 W). Consult your facilities department regarding ventilation requirements for this level of heat output.
Emission/ Immunity Statement	For our European customers, any product marked with the CE label meets the European requirements for emission and immunity as defined in the EMC Directive 89/336/EEC. This product has been evaluated to the Standard for Emissions for Industrial Scientific Equipment (EN 55011 – Class A), and to the Standard for Generic Immunity (EN50082-1).

Laboratory Ventilation

Venting This instrument produces toxic gaseous waste that must be ventilated either through a fume hood to a duct or directly to the duct. The fume hood or duct must be located within 5 m (15 ft.) of the instrument. The information presented here reflects U.S. regulations and practices for venting waste from Applied Biosystems instruments to a fume hood or to a duct.

! WARNING ! Some Applied Biosystems instruments use chemicals that are hazardous. Always mix and prepare hazardous materials under an operating fume hood.

CAUTION Dispose of all waste in accordance with all applicable local, state, and federal environmental health and safety regulations and laws.

Fume Hood Following are important points about the fume hood:

- ◆ The fume hood should operate continuously, including nights and weekends, because vented waste bottle contents can escape to surroundings.
 - ◆ The fume hood must be constructed of materials that are compatible with the waste materials/chemicals being generated or exhausted.
 - ◆ The fume hood should be located away from air currents generated by air conditioning ducts, fans, windows, doors, and moving equipment and persons.
 - ◆ The fume hood exhaust vent should be located where gaseous waste cannot be drawn back in the building.
 - ◆ A sign or label should be present that shows where to locate the fume hood sash so as to give an average flow of 100 ft./min. (linear) face velocity. The minimum velocity flow at any point in the hood must be 80 ft./min. (linear). The maximum flow must not exceed 125 ft./min. (linear).
 - ◆ The fume hood must meet local, state, and federal health and safety requirements. Refer to current fume hood standards established by the American Society of Heating, Refrigeration, Air conditioning engineers (ASHRAE), American Conference of Governmental Industrial Hygienists (ACGIH), and Occupational Safety and Health Agency (OSHA).
 - ◆ Check and record face velocity at least yearly.
 - ◆ Inspect and maintain exhaust system, including fans and motors at least yearly.
-

Duct System It is important that the duct system:

- ◆ operate continuously, including nights and weekends, because vented waste bottle contents can escape to surroundings.
- ◆ be constructed of PVDF tubing or other materials compatible with the waste materials being generated.
- ◆ not come in contact with strong oxidizers, bases, or other chemicals that are incompatible with gaseous waste.
- ◆ allow vapor or gas movement of 1000–2000 ft./min. (5.5–11 m/s).

! WARNING ! Do not connect the waste vent to a ductless hood or to a system that purifies, filters air, and returns it to the room.

Instrument Safety

2

Instrument Safety

Safe Operation Before operating the instrument, read the information in this section concerning hazards and potential hazards. Ensure that anyone involved with the operation of the instrument is instructed in both general safety practices for laboratories and specific safety practices for the instrument.

User Attention Words Four user attention words appear in the text of all documentation for Applied Biosystems products. Each word implies a particular level of observation or action that is significant to user safety or proper instrument operation.

Note Calls attention to useful information.

IMPORTANT Indicates information that is necessary for proper instrument operation.

CAUTION Cautions the user that a potentially hazardous situation could occur, causing injury to the user or damage to the instrument, if this information is ignored.

! WARNING ! Warns the user that serious physical injury or death to the user or other persons could result if these precautions are not taken.

Instrument Labeling

Safety Labels Safety labels are located on the instrument. Each label consists of a Signal Word panel and a Message Panel. A Safety Alert Symbol indicates a potential personal safety hazard. If multiple hazards exist, the signal word corresponding to the greatest hazard is used.

Signal Words

- ◆ **DANGER** indicates an imminently hazardous situation that, if not avoided, will result in death or serious injury (most extreme).
- ◆ **CAUTION** indicates a potentially hazardous situation that could result in minor or moderate injury to the user or damage to the instrument.
- ◆ **WARNING** indicates a potentially hazardous situation that could result in death or serious injury.

Labels That May Be Found On Your Instrument

The following Danger, Caution, and Warning labels, listed in English and French below, may be found on your instrument.

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Safety Alert Symbols

Electrical Symbols The following chart is an illustrated glossary of all electrical symbols that are used on Applied Biosystems instruments. Whenever such symbols appear on instruments, please observe appropriate safety procedures.

	This symbol indicates the on position of the main power switch.
	This symbol indicates the off position of the main power switch.
	This symbol indicates the on/off position of a push-push main power switch.
	This symbol indicates that a terminal may be connected to another instrument's signal ground reference. This is not a protected ground terminal.
	This symbol indicates that this is a protective grounding terminal that must be connected to earth ground before any other electrical connections are made to the instrument.
	A terminal marked with this symbol either receives or delivers alternating current or voltage.
	A terminal marked with this symbol can receive or supply an alternating and a direct current or voltage.
	CAUTION This symbol indicates the presence of high voltage and warns the user to proceed with caution.
	CAUTION This symbol alerts you to consult the manual for further information and to proceed with caution.

Non-electrical Symbols The following is an illustrated glossary of all non-electrical safety alert symbols found on Applied Biosystems instruments.

	CAUTION This symbol illustrates a heater hazard. Proceed with caution when working around these areas to avoid being burned by hot components.
	This symbol indicates that a laser is present inside the instrument.

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Input/Output Connections

Location The exact locations of the input/output connections on the PROCISE/PROCISE cLC Protein Sequencer are shown in Figure 2-1. The input/output connections are labeled on the instrument.

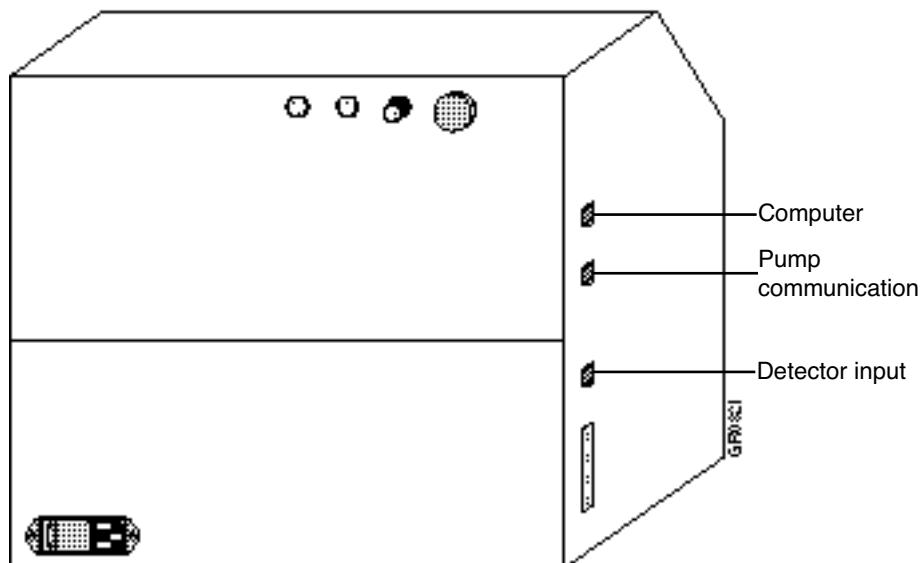
Figure 2-1 Input/Output connections on the PROCISE/PROCISE cLC Protein Sequencer

Gas Safety

Pressurized Gas Safety **! WARNING !** Pressurized gas cylinders are potentially explosive if not handled properly, and can cause great damage and severe injuries.

This instrument uses a pressurized gas cylinder for operation. Some kinds of rapidly leaking gas can displace normal atmosphere and cause suffocation. If knocked over, a pressurized gas cylinder can explode. Keep in mind the following when working with gas cylinders:

- ◆ Screw the gas cap on tightly when not in use and when transporting.
- ◆ Attach the cylinder firmly to a wall or gas cylinder cart by means of an approved strap or clamp.

Gaseous Waste Safety **! WARNING !** Handle all liquid, solid, and gaseous waste as potentially hazardous. Gaseous waste can cause injury, illness, or death.

Instruments that use gas cylinders for operation emit gaseous waste. Follow carefully the instructions provided in this chapter in the section entitled, "Disposing of Gaseous Waste and Exhaust," and read the MSDSs and Waste Profiles in Chapter 3, "Chemical Safety."

Keep in mind the following when working with instruments that use gas cylinders:

- ◆ Some gaseous waste is toxic.
 - ◆ Vent according to instructions in the guide and all appropriate laws and regulations.
-

Disposing of Gaseous Waste

Instrument Waste and Exhaust System The waste and exhaust system of the instrument is composed of a 3/8 in.-o.d. (0.95 cm) exhaust line attached to a 2.5-gal. (9.5 L) polyethylene bottle in a secondary container. This exhaust line conducts the waste fumes and gasses that enter the waste manifold to a fume hood or duct system for disposal.

Note In the 432 Peptide Synthesizer, the waste and exhaust system has a 1/4 in.-o.d. exhaust line attached to a 1-gal (4 L) polyethylene bottle.

CAUTION Be sure to place the waste bottle in a secondary container to minimize the possibility of leaks.

The fume hood or duct system must be operating whenever the instrument power is on or when there is waste in the waste container.

Exhaust Line Tubing

- ◆ Exhaust line tubing should be polypropylene tubing of the shortest possible length and the straightest possible run. Tubing length should not exceed 15 ft (4.5 m).
- ◆ The exhaust line tubing should not have low points that can trap residue or condensation.
- ◆ The exhaust line tubing should be fastened securely. Use fasteners of polypropylene or Teflon. Do not use brass as it corrodes. Be careful not to puncture tubing.
- ◆ The tubing should be located away from sources of potential damage such as contact, heat, or flame.
- ◆ The tubing end should be placed as far as possible into the duct, canopy, or hood.
- ◆ The open end of the tubing should not face into oncoming air movement through the duct or canopy.

Connecting the Exhaust line Connect the exhaust line from the instrument to the waste bottle so that it drops vertically. This prevents liquid and waste from accumulating and blocking the flow.

Make sure that the line leads to the fume hood or duct in a straight and upward direction. The line must not create low points by dipping in a downward direction because condensation may collect and prevent proper flow through the instrument.

Pressurized Gas and Accessories Needed but Not Supplied

Gas Cylinder Each laboratory must supply the required gas cylinder and accessories for installation. This instrument requires one size 1-A argon gas cylinder that holds approximately (7.2 m³) 257 ft.³ of gas when full.

Use only prepurified argon of 99.998% or greater purity. An additional argon cylinder should be ordered in advance for the ongoing operation of the instrument.

CAUTION Damage to the instrument and its products can result from using impure argon, gasses other than argon, or an inadequate amount of argon.

Pressure Regulator A two-stage pressure regulator with dual gauges (output range: 0–80 psi, high pressure range: 0–3000 psi), and a Compressed Gas Association (CGA) 580 cylinder adapter with a needle-type shutoff valve on the exit side is required. The needle valves should have Swagelok- type end-fittings ready for connection to .25-in. (6.35 mm) o.d. tubing.

CAUTION Do not allow the cylinder pressure to drop below 300 psi. Chemical vapors and/or liquids can backflow into the pressure regulator resulting in serious damage to the instrument.

The second stage output of the regulator should be set at approximately 15 psi (103.5 kPa).

Attaching the Cylinder Attach the pressurized gas cylinder firmly to a wall or gas cylinder cart by means of approved straps or clamps.

! WARNING ! Pressurized gas cylinders are explosive. Attach the pressurized gas cylinder firmly to a wall or bench by means of approved straps or clamps. Always cap the gas cylinder when not in use.

Chemical Safety

3

Introduction

Please read... This chapter contains detailed information for this instrument about:

- ◆ hazardous chemicals used (if any)
- ◆ hazardous waste produced and how to handle

Where appropriate, detailed Waste Profiles and Material Safety Data Sheets (MSDSs) are provided to ensure correct and safe operation of the instrument.

Applied Biosystems assumes that all operations in your laboratory will be conducted in accordance with safety practices detailed in Waste Profiles and MSDSs and with any applicable laws and regulations.

Hazardous Chemicals

Overview This instrument uses chemicals that are hazardous, chlorinated, organic liquids.

! WARNING ! Hazardous chemicals can cause injury, illness, or death. Handle all chemicals as potentially hazardous.

Chemicals are classified as hazardous when they are physically hazardous or if they may cause health hazards upon acute or chronic exposure.

- ◆ Physically hazardous chemicals are materials that are flammable, combustible, compressed gasses, explosives, oxidizers, organic peroxides, pyrophoric, reactive or unstable, or water reactive.
- ◆ Chemicals that may cause health hazards are materials that are carcinogens, toxic or highly toxic, reproductive toxins, irritants, corrosives, sensitizers, hepatogens, nephrotoxins, agents which act on the hematopoietic systems, and agents which damage the lungs, skin eyes, or mucous membranes.

Handling Hazardous Chemicals Important requirements for handling hazardous chemicals include:

- ◆ Read all applicable MSDSs before handling hazardous chemicals.
 - ◆ When replacing reagents, always install the new bottle on the instrument. Do not top off old bottles. Some chemicals reduce the integrity of glass bottles and repeated use beyond six weeks may result in the bottle fracturing when it is pressurized during operation.
 - ◆ Provide secondary containment for all reagent bottles.
 - ◆ Do not store chemicals in direct sunlight or heat (on or off the instrument).
-

Hazardous Waste

Overview The PROCISE/PROCISE cLC Protein Sequencer generates hazardous, chlorinated, organic liquid and gaseous waste.

! WARNING ! CHEMICAL WASTE HAZARD. Waste produced by this instrument is potentially hazardous and can cause physical injury, illness, or death. Dispose of the contents of the waste tray and waste bottle in accordance with all applicable health and environmental laws and regulations.

Handling Hazardous Waste

The following are important requirements for handling hazardous waste:

- ◆ Ensure that the waste container is correctly installed.
 - ◆ Always handle hazardous materials beneath a fume hood that is connected in accordance with all installation requirements.
 - ◆ During transporting, ensure that waste containers are tightly sealed with the cap provided.
 - ◆ Read the Instrument Waste Profiles in this chapter before handling or disposing of hazardous waste.
 - ◆ Read all applicable Material Safety Data Sheets before handling or disposing of hazardous waste.
 - ◆ Dispose of hazardous waste in accordance with all applicable laws and regulations.
-

Storing Hazardous Waste

! WARNING ! Never collect or store waste in a glass container because of the risk of breaking or shattering.

The following are guidelines for storing hazardous waste:

- ◆ Do not store waste for long periods of time.
 - ◆ Do not store large amounts of waste in the lab.
 - ◆ Do not store waste in direct sunlight or heat (on or off instrument).
-

Instrument Waste Profile Overview

General Information

The Instrument Waste Profile provides essential information about hazardous waste. Read the Instrument Waste Profile and all applicable MSDSs before handling or disposing of waste. The Instrument Waste Profile is NOT a substitute for MSDS information. See specific MSDS sheets in this chapter for chemical constituent, health, toxicological, and hazard information. For a complete list of acronyms and abbreviations used in Waste Profiles, see Appendix A.

Waste Profile Sections

Information in the Instrument Waste Profiles is divided into ten sections. The names of the ten sections are listed below.

Topic	Section
Identification	1
Approximate Composition	2
Physical Data	3
Fire and Explosion Hazard Data	4
Health Hazard Data	5
Reactivity Data	6
Spill or Leak Procedures	7
Special Protective Equipment	8
Special Precautions	9
Additional Information	10

PROCISE/PROCISE cLC Protein Sequencing System Waste Profile

Emergency Phone Numbers:(USA) 650-570-6667 (UK) 0925-825650

section 1 - IDENTIFICATION

The liquid waste from PROCISE cLC is collected in a one liter rectangular glass waste bottle located behind the right-hand door of the instrument at the lower right. A spring-loaded aluminum cap holds the bottle in place. The cap is simply lifted up to remove the waste bottle. This waste is a complex mixture of reagents and solvents which may have properties of greater hazard than the individual components by themselves. The waste consists of two phases. The top non-aqueous phase is highly flammable. The bottom aqueous phase contains residual trifluoroacetic acid and is corrosive.

HANDLE THIS MATERIAL WITH EXTREME CAUTION! DO NOT DISPOSE OF THIS WASTE IN SINKS OR DRAINS! THIS MATERIAL SHOULD BE DISPOSED OF AS A REGULATED HAZARDOUS WASTE!

Section 2 - APPROXIMATE COMPOSITION

PROCISE

Material	%	TLV	PEL	STEL	CAS#
Ethyl Acetate (EtAc)	40	400 ppm	400 ppm		141-78-6
Butyl Chloride	31	N/A	N/A		N/A
Water	17				
Acetonitrile	4	40 ppm	40 ppm	60 ppm	75-05-8
Heptane	3	400 ppm	400 ppm	500 ppm	142-84-5
Trifluoroacetic acid (TFA)	2	N/A	N/A		N/A
Phenylisothiocyanate (PITC)	1	N/A	N/A		N/A
N-methylpiperidine	1	N/A	N/A		626-67-5
Methanol	1	200 ppm	N/A	250 ppm	67-56-1

N/A = Not Available

NFPA Combined Health Rating is: 3 (Health); 4(Flammability); 3 (Reactivity)

PROCISE cLC

Material	%	TLV	PEL	STEL	CAS#
Butyl Chloride	45	N/A	N/A		N/A
Ethyl Acetate (EtAc)	26	400 ppm	400 ppm		141-78-6
Water	16				
Heptane	6	400 ppm	400 ppm	500 ppm	142-84-5
Acetonitrile	4	40 ppm	40 ppm	60 ppm	75-05-8
Trifluoroacetic acid (TFA)	2	N/A	N/A		N/A

PROCISE CLC

Phenylisothiocyanate (PITC)	1	N/A	N/A		N/A
N-methylpiperidine	1	N/A	N/A		626-67-5
Methanol	1	200 ppm	N/A	250 ppm	67-56-1

N/A = Not Available

NFPA Combined Health Rating is: 3 (Health); 4(Flammability); 3 (Reactivity)

Section 3 - PHYSICAL DATA

BOILING POINT	760 mm: N/A
SPECIFIC GRAVITY	(H ₂ O = 1): 1.01
VOLATILITY	(vol%): 98-99%
APPEARANCE AND ODOR:	Colorless two-phase liquid with a fruity, vinegar odor.
FREEZING POINT	N/A
pH RANGE	0-2 (bottom phase)
SOLUBILITY IN H ₂ O	Top 5%; btm 100%

Section 4 - fire and explosion hazard data (Acetonitrile data only!)

FLASH POINT (Closed Cup):	(Ethyl Acetate only) -5 °C (23 °F)
FLAMMABLE LIMITS:	(Ethyl Acetate only) 2.2% LEL; 11% UEL
FIRE EXTINGUISHING MEDIA:	Dry chemical, alcohol foam, carbon dioxide or Halon
SPECIAL FIRE FIGHTING PROCEDURES:	Use self-contained breathing apparatus and protective clothing to prevent skin and eye contact.

Section 5 - HEALTH HAZARD DATA

EXPOSURE LIMITS:	See Section II.
	For Ethyl Acetate, the IDLH level is 10,000 ppm. For Heptane, NIOSH recommends a TWA of 85 ppm and a ceiling level of 440 ppm. The STEL is 500 ppm and the IDLH level is 4250 ppm. For Acetonitrile, the STEL is 60 ppm and the IDLH level is 4,000.
EFFECTS OF ACUTE OVEREXPOSURE:	
SWALLOWING	Harmful if swallowed! Causes severe irritation of eyes, nose, and throat, Higher concentrations may cause liver and kidney damage, unconsciousness, and death.
SKIN	May cause severe irritation or burns. Allergic skin sensitization may also occur.

EXPOSURE LIMITS:	See Section II.
INHALATION	May cause irritation of eyes, nose, throat, and lungs. Higher concentrations may cause pulmonary edema, unconsciousness, and death.
EYES	May cause severe irritation or burns.
EMERGENCY AND FIRST-AID PROCEDURES	If conscious, give large quantities of water immediately and induce vomiting. Get medical attention immediately. Do not induce vomiting if unconscious.
SKIN:	Remove contaminated clothing. Flush the contaminated area with water and wash with mild soap or detergent. Get medical attention.
INHALATION:	Provide fresh air and rest. If breathing is difficult, provide oxygen and get medical attention immediately.
EYES:	Flush eyes immediately with large amounts of water for at least 15 minutes. Get medical attention.

Section 6 - REACTIVITY DATA

STABILITY	Stable
INCOMPATIBILITY:	Contact with strong oxidizing agents or concentrated acids or bases may cause fire or explosion.
HAZARDOUS COMBUSTION OR DECOMPOSITION PRODUCTS:	Burning may release toxic vapors and gases including phosgene, hydrogen chloride, hydrogen fluoride, carbon monoxide, and oxides of nitrogen.
HAZARDOUS POLYMERIZATION:	Will not occur.

Section 7 - SPILL OR LEAK PROCEDURES

STEPS TO BE TAKEN:	Avoid inhalation and skin contact. Wear protective clothing. Ventilate area of spill or leak. Remove all ignition sources. Small quantities may be collected with absorbent towels or pads and removed to a well-ventilated area away from ignition sources. Larger amounts (1 liter or more) may be collected with an inert absorbent (kitty litter or similar material) or commercially available spill pillows designed for solvent collection. This waste material must not be allowed to enter confined spaces (such as a sewer) because of the possibility of an explosion.
WASTE DISPOSAL:	This instrument waste solution should be disposed of as a regulated hazardous waste by a properly-permitted hazardous waste management facility in accordance with federal, state and local regulations. Recommend disposal methods include high temperature incineration and solidification for secure chemical landfill disposal.

Section 8 - SPECIAL PROTECTIVE EQUIPMENT

RESPIRATORY PROTECTION:	An MSHA- or NIOSH-approved respirator for organic vapors is recommended. A supplied-air or SCBA respirator is recommended for high vapor concentration and emergency situations.
VENTILATION:	Handle within a well-ventilated area. Minimize open exposure to air.
PROTECTIVE GLOVES:.	Butyl rubber or similar type resistant to chlorinated solvents.
EYE PROTECTION:	Safety glasses with side shields, monogoggles, or face shield.
OTHER PROTECTIVE EQUIPMENT:	As necessary to prevent skin contact.

Section 9 - SPECIAL PRECAUTIONS

PRECAUTIONS TO BE TAKEN:	Handle as a flammable, poisonous liquid. Maintain adequate ventilation at all times. Do not breathe vapor. Do not get in eyes, on skin, or on clothing. Accidental contact should be washed immediately. Keep away from heat, sparks, and flame. Spill collection materials, eye wash, and safety shower should be in area of use.
OTHER:	This waste solution has strong solvent properties and will attack many forms of rubber, plastics, coating, and finishes.

Section 10 - ADDITIONAL INFORMATION

When not directly attached to the instrument, this waste material should be stored in a secure, well-ventilated location suitable for flammable materials. Store away from light, heat, or potential ignition sources. Contact the appropriate state hazardous waste regulatory agency for proper disposal procedures and lists of registered service companies.

THIS WASTE MATERIAL IS HAZARDOUS AND SHOULD ONLY BE HANDLED BY PERSONS THOROUGHLY TRAINED IN HAZARDOUS MATERIALS HANDLING PROCEDURES!

Material Safety Data Sheets

Overview Material Safety Data Sheets provide information about physical characteristics, health hazards, safety precautions, first aid, spill cleanup and disposal procedures. Read the MSDSs before handling reagents or interacting with the instrument. For a complete list of acronyms and abbreviations used in MSDSs, see Appendix A.

Contents of MSDSs Each MSDS is divided into 16 sections, although not all sections apply to every chemical. The titles of the 16 sections are listed in the table below.

Title / Topic	Section #
Chemical Product and Company Identification	1
Composition, Information on Ingredients	2
Hazards Identification	3
First Aid Measures	4
Fire Fighting Measures	5
Accidental Release Measures	6
Handling and Storage	7
Exposure Controls, Personal Protection	8
Physical and Chemical Properties	9
Stability and Reactivity	10
Toxicological Information	11
Ecological Information	12
Disposal Considerations	13
Transport Information	14
Regulatory Information	16
Other Information	12

How to Get Extra Copies of MSDSs Additional copies of MSDSs for chemicals manufactured by Applied Biosystems are available at no extra cost. Call your local representative or access the Applied Biosystems World Wide Web site at www.appliedbiosystems.com/techsupport

MSDSs for the PROCISE/PROCISE cLC Protein Sequencer

MSDSs for Applied Biosystems Chemicals

The MSDSs for chemicals manufactured by Applied Biosystems and used with the PROCISE/PROCISE cLC Protein Sequencer follow this section. They are:

- ◆ Biobrene Plus (P/N 902420)
- ◆ Methanol (P/N 901834)
- ◆ Premix buffer concentrate (P/N 902173)
- ◆ R1, 5% PITC in heptane (P/N 901824)
- ◆ R2, N-Methylpiperidine/water/MeOH (P/N 902533)
- ◆ R3, TFA (P/N 901812)
- ◆ R4A, 25% TFA in water (P/N 901816)
- ◆ R5, acetonitrile (P/N 901817)
- ◆ S1, n-Heptane (P/N 901818)
- ◆ S2B, ethylacetate (P/N 901813)
- ◆ S3, n-Butyl Chloride (P/N 901814)
- ◆ S4B, 20% acetonitrile in water (P/N 901829)
- ◆ Solvent A3 (3.5% THF in water) (P/N 902209)
- ◆ Solvent B2 (acetonitrile/2-Propanol) (P/N 902675)
- ◆ S4C, 10% Acetonitrile in water (P/N 901829)

MSDSs for Chemicals Not Supplied by Applied Biosystems

For chemicals required for this instrument but not manufactured or sold by Applied Biosystems, obtain the MSDSs from their manufacturers. Brief warnings for these chemicals are listed below.

! WARNING ! CHEMICAL HAZARD. Acetone is a flammable chemical, irritating to respiratory system if inhaled. It can cause liver , kidney, and eye damage. Wear protective clothing, chemical resistant gloves, and protective eyewear.

! WARNING ! CHEMICAL HAZARD. Methanol is a flammable chemical, irritating to respiratory system if inhaled. It can cause liver , kidney, and eye damage. Wear protective clothing, chemical resistant gloves, and protective eyewear.

! WARNING ! CHEMICAL HAZARD. Sodium Dodecyl Sulfate (SDS) is a toxic chemical that is harmful to lungs and internal organs if swallowed. Contact with eyes may cause serious damage. Wear protective clothing, chemical-resistant gloves, and protective eyewear.

Acronyms and Abbreviations

Acronyms and Abbreviations Used in Waste Profiles and MSDSs

Introduction Waste Profiles and MSDSs use acronyms and abbreviations for certain organizations, government regulations, common scientific terminology, units of measurement, and chemicals.

The following tables contain lists of acronyms and abbreviations commonly used in Waste Profiles and MSDSs.

Organizations, Regulations, and Scientific Terminology

Acronyms and abbreviations for organizations, government regulations, and scientific terminology.

Table 1-1 Acronyms and abbreviations

Term	Explanation
ACGIH	American Conference of Governmental Industrial Hygienists
CAS#	Chemical Abstract Service Reference Number for Specific Pure Chemical
cc	Closed cup testing of flash point
CFR	Code of Federal Regulations. Regulations published by the United States Government
CERCLA	Comprehensive Environmental Response, Compensation and Liability Act (Superfund) is a federal law administered by EPA
DFG MAK	Federal Republic of Germany's Maximum Contamination Value in the workplace (similar to PEL in the USA)
DOT	United States Department of Transportation, regulates transportation of hazardous material (USA)
EPA	United States Environmental Protection Agency, regulates use, disposal, or emission of hazardous material (USA)
IDLH	Immediate Danger to Life and Health
LC _{LO}	Lowest published lethal concentration
LC ₅₀	Lethal concentration in air that kills 50% of a specified population
LD ₅₀	Lethal dose that kills 50% of a specified population
LEL	Lower Explosion Limit
MSHA	Mine Safety and Health Administration, recommends respirators
NFPA	National Fire Protection Association, publishes recommended regulations for local or state governments in the United States (Hazardous rating system developed by this Association)

Table 1-1 Acronyms and abbreviations *(continued)*

Term	Explanation
NIOSH	National Institute of Occupational Safety and Health (USA) recommends exposure levels and respirators
oc	Open cup testing for flash point
OSHA	Occupational Safety and Health Administration (USA), sets chemical exposure levels.
PEL	Permissible Exposure limit. The federal OSHA limit, usually expressed as TWA for an 8-hour work shift.
PPM	Parts Per Million
Prop 65	A California Law requiring warnings for chemicals that are known to the state to be carcinogenic or cause reproductive harm.
RCRA	Resource Conservation and Recovery Act
RTECS	Registry of Toxic Effects of Chemical Substances
SARA	Superfund Amendments and Reauthorization Act, a federal act administered by EPA.
SCBA	Self-Contained Breathing Apparatus
STCC	Standard Transportation Commodity Code
STEL	Short Term Exposure Level, published by ACGIH
TC _{LO}	Lowest published toxic concentration
TLV	Threshold Limit Value. The ACGIH-recommended TWA, usually for an 8-hour work shift
TWA	Time Weighted Average
UEL	Upper Explosive Limit
u or U	Unknown
UN	United Nations. This designation identifies hazardous chemicals in the process of world-wide transportation.

Units of Measurement

The following are abbreviations for units of measurement.

Table 1-2 Abbreviations for units of measurement

Abbrev.	Unit of Measurement
#	number
°C	degrees Celsius
°F	degrees Fahrenheit
μL	microliter
μm	micron
μmol	micromole
AUFS	absorbency units full-scale
ft.	foot
i.d.	inside diameter
in.	inch
L	liter
m	meter

Table 1-2 Abbreviations for units of measurement

Abbrev.	Unit of Measurement
mg	milligram
mL	milliliter
mm	millimeter
o.d.	outside diameter
P/N	Applied Biosystems part number
psi	pounds per square inch
sec	second
V	volt

Chemicals The following are abbreviations for chemicals.

Table 1-3 Abbreviations for chemicals

Abbrev.	Definition
A	adenine
AA	amino acid
1Ac	acetyl
AcI	acetylimidazole
Acm	acetamidomethyl
Ac ₂ O	acetic anhydride
ACN	acetonitrile
ACT	activator vessel
BHA resin	benzhydramine resin
t-Boc	tert-butyloxycarbonyl
Bzl	benzyl
Br-Z	2-bromobenzyloxycarbonyl
t-Bu	tert-butyl
C	cytosine
CHO	formyl
CH ₃ Bzl	4-methylbenzyl
CH ₃ O Bzl	4-methoxybenzyl
Cl-Z	2-chlorobenzyloxycarbonyl
CPG	Controlled Pore Glass
DCA	dichloroacetic acid
DCC	dicyclohexylcarbodiimide
DCM	dichloromethane
DCU	dicyclohexylurea
DIEA	diisopropylethylamine
DMAP	4-dimethylaminopyridine
DMF	dimethylformamide
DMSO	dimethylsulfoxide
DNA	deoxyribonucleic acid

Table 1-3 Abbreviations for chemicals *(continued)*

Abbrev.	Definition
Dnp	2,4-dinitrophenyl
Et	ethyl
EtOH	ethanol
Fmoc	9-fluorenylmethyloxycarbonyl
G	guanine
HBTU	2-(1 H-benzotriazol-1-yl)-1,1,3,3-tetramethyl-uronium hexafluorophosphate
HLP	high loaded polystyrene
HMP resin	p-hydroxymethylphenoxymethyl-polystyrene resin
HOAc	acetic acid
mBHA resin	4-methylbenzhydramine resin
MeOH	methanol
Mob	4-methoxybenzyl
Mtr	4-methoxy-2,3,6-trimethyl-benzene sulfonyl
Mts	mesitylene-2-sulfonyl
NMI	1-Methylimidazole
NMP	N-Methylpyrrolidone, N-methyl-2-pyrrolidone
OBt	ethyl ester
OMe	methyl ester
PAM resin	phenylacetamidomethyl resin
PEG	polyethylene glycol
RV	reaction vessel
SSPS	solid-phase peptide synthesis
T	thymine
TETD	tetraethylthiuram disulfide
TFA	trifluoroacetic acid
TFMSA	trifluoromethane sulfonic acid
THF	tetrahydrofuran
Tos	4-toluenesulfonyl (tosyl)
Tri	trityl
U	uracil
Z	benzyloxycarbonyl

Headquarters

850 Lincoln Centre Drive
Foster City, CA 94404 USA
Phone: +1 650.638.5800
Toll Free: +1 800.345.5224
Fax: +1 650.638.5884

Worldwide Sales Offices

Applied Biosystems vast distribution and service network, composed of highly trained support and applications personnel, reaches into 150 countries on six continents. For international office locations, please call our local office or refer to our web site at www.appliedbiosystems.com.

www.appliedbiosystems.com

Applera Corporation is committed to providing the world's leading technology and information for life scientists. Applera Corporation consists of the Applied Biosystems and Celera Genomics businesses.

Printed in the USA, 01/2002
Part Number 904203E

an **Applera** business